



**UNIVERSITY OF
BIRMINGHAM**

**Year 3 Laboratory:
Interactions of Neutrons with Matter**

School of Physics and Astronomy
University of Birmingham

Josh Wainwright
UID:1079596

Demonstrators: Dr Peter Jones, Dr Garry Tungate,
Dr Carl Wheldon, Dr Neil Curtis
Date: December 2012

Contents

1	Introduction	2
2	Neutron Binding Energy	2
2.1	Sodium Iodide Detector	2
2.2	Calibration	3
2.3	Energy Resolution	4
2.4	Neutron Binding Energy	5
3	Neutron Attenuation In Water	6
3.1	Boron Tri-Flouride Detector	6
3.1.1	Wall Effect	7

1 Introduction

This experiment involves examining the interactions of neutrons with matter through two experiments using two different detectors. The first detector, a sodium iodide scintillation detector will be used to measure the binding energy of the neutron, and secondly, using a boron-trifluoride counter, measure the moderating properties of water on fast neutrons as a function of radius from a central source.

2 Neutron Binding Energy

2.1 Sodium Iodide Detector

The sodium iodide detector is a scintillation detector that uses a NaI crystal doped with thallium to create scintillation photons when a neutron hits the crystal. The thallium is added as an activator to provide extra available energy levels for the photons in the crystal to occupy. In the pure state, the NaI has a large forbidden band gap between the conduction and valence bands, shown in figure 1. The activator reduces this gap so that an emitted photon is in the range that the photomultiplier tube is sensitive to, typically visible.

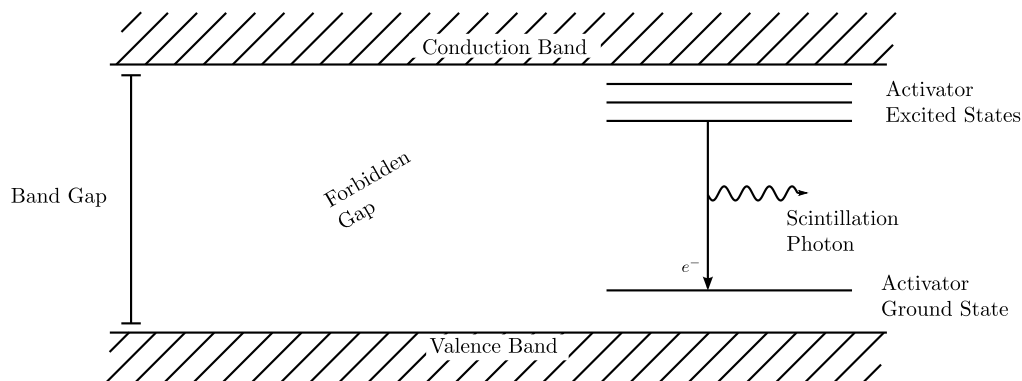


Figure 1: The element thallium is added to the pure crystal of sodium iodide to provide extra energy levels for the electrons, which sit inside the forbidden gap of the crystal lattice structure. These allow a photon of a useful wavelength to be emitted when electron decays to a lower energy.

When a gamma photon enters the detector it ionises a molecule of NaI creating excited states in the crystal that then decay via visible photons. The number of photons that are produced is proportional to the energy of the incident photon, as the energy of the gamma ray increases, more scintillation photons are produced. Thus, the energy of the original can be found by counting the number of photons produced.

These scintillation photons are used to create electrons which are more easily counted. At the photo-cathode, the photons release electrons via the photoelectric effect. The electrons are then passed through the photomultiplier tube which increases the number of photons linearly by accelerating them through a high potential difference so they can knock more electrons out from the dynodes (see figure 2). This results in a stream of electrons that can be measured, and which is proportional to the original energy of the photon via the number of scintillation photons, the number of electrons and the multiplication factor of the increase of electrons from the PM tube. Since each of these steps is a linear relationship, a calibration setting can be taken using known results, and a simple linear equation relating the reading taken to the energy of the photon found.

The detector requires a high voltage which is placed across the dynodes and the anode to accelerate the electrons once they have entered the PM tube. Using this method of acceleration


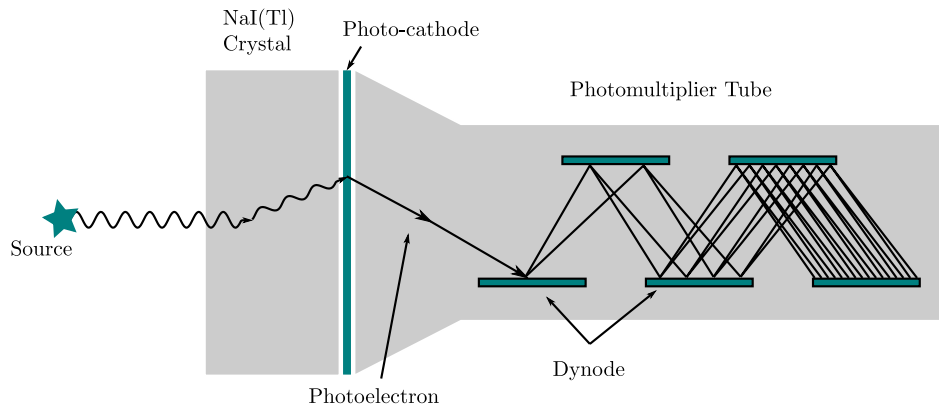


Figure 2: The sodium iodide detector is a commonly used gamma ray detector that uses the relationship between the energy of a photon and the ionisation power that it has to measure its energy via a calibration equation.

and re-acceleration, the photomultiplier tube is able to increase the number of electrons by a factor of roughly 10^5 depending on the voltage bias.

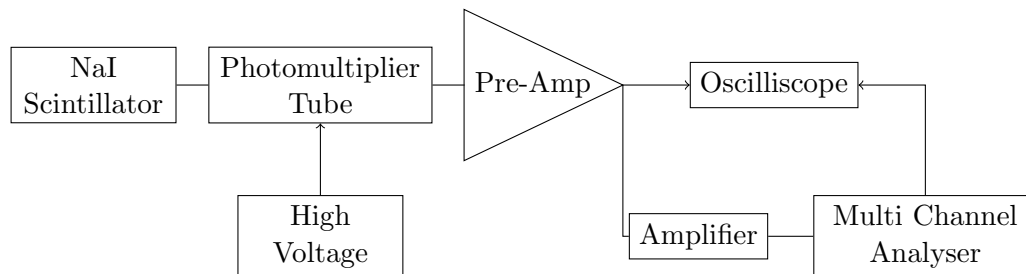

2.2 Calibration

Before any data can be read from the spectra taken from the detector, the data has to be calibrated. To do this, a set of known samples with well defined peaks is used to get a relation between the channel number of an observed peak and its energy. Once this is done, a relation can be found that relates the channel number to energy for use on later measurements.

The elements that were used included cobalt, barium, caesium, europium and americium. Each of these elements were available to measure, and each has recorded and accepted peaks that are at well defined positions. The data was taken from the National Nuclear Data Centre (NNDC). The table below, table 1, shows the data that we collected from these sources. The channel number refers to the position in the spectrum that the peak was observed and the expected peak is the energy of the peak that the respective peak is estimated to correspond to.

From these, the graph in figure 3 can be drawn which shows the relation that can be used as a calibration between energy and channel number. The exact value of this relation, intercept and gradient, are subject to change since the calibration is retaken at the start of each session to account for changes in the experimental set-up and the radioactivity of the source which will decrease with time.

An example of the spectra from the sources is shown in figure 4. This shows the spectrum of ^{60}Co which has been fitted with the Gaussian peaks to find the position and FWHM using the “root” package.

Element	Channel Number	Expected Peak Energy (keV)
Cobolt ^{60}Co	369.4	1332.0
	324.9	1173.0
Barium ^{133}Ba	102.9	356.0
	86.08	302.9
Caesium ^{137}Cs	188.4	661.7
Europium ^{152}Eu	99.84	344.3
	222.0	778.9
	274.9	964.1
	311.0	1121/1086
	394.1	1408.0

Table 1: Several elements were used to measure a calibration relation that links the observed channel number to the energy.

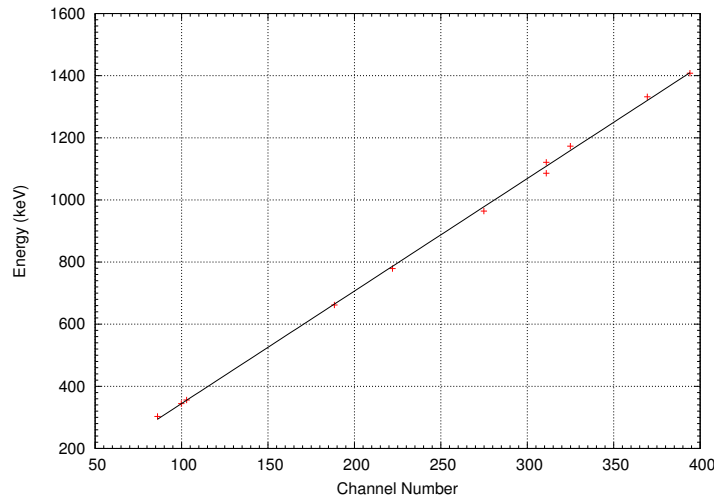


Figure 3: The sodium iodide detector provides a channel number readout that is directly proportional to the energy of the corresponding reaction event. Thus, using known data points, a relation between channel number and energy can be found.

2.3 Energy Resolution

Scintillation detectors, such as the NaI detector in use here, have a relatively poor energy resolution. Statistical broadening of the peaks is usually the dominant form of broadening. When a reaction takes place, and so energy is deposited, in the scintillator crystal, the efficiency can be as low as around 10%. This means that a far smaller amount of energy is actually transferred to the crystal than the energy of the incoming particle.

When the reaction takes place in the crystal, a number of photons, proportional to the energy deposited, are produced. But some of them are lost before they can reach the photomultiplier tube further reducing the measured number of photons. This also means that there is a minimum in the number of photons just before the photomultiplier tube. Since the number of arriving photons will be subject to statistical fluctuations, the standard deviation in the photons reaching the PM tube is given by equation 2.1,

$$\text{Standard Deviation, } \sigma = \sqrt{\text{Number of photons reaching the PM tube}} \quad (2.1)$$

Since this is proportional to the energy of the original photon that caused the reaction, it can

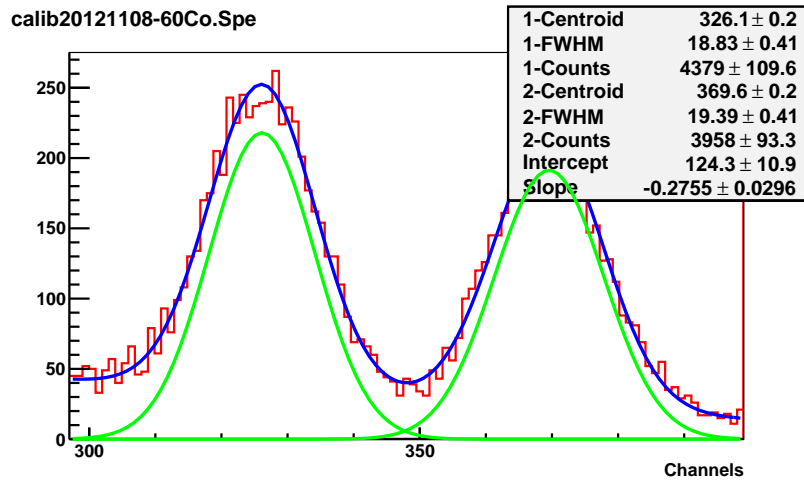


Figure 4: An example of the peak for the radioactive element cobalt-60. These two peaks are well defined and have well known values. This means that they are good for using as calibration points for the detector.

be shown that

$$R = \frac{\text{FWHM}}{H_0} \quad (2.2)$$

$$= k \frac{\sqrt{E}}{E} = \frac{k}{\sqrt{E}} \quad (2.3)$$

Where FWHM is the full width at half maximum of any given peak, H_0 is the height of the peak and k is a constant of proportionality. This means that the resolution can be easily observed by plotting a graph of $\log k$ against $\log E$, which should give a slope of $-\frac{1}{2}$ in an ideal case.

The resolution of the NaI detector is use was measured using the method described above.

2.4 Neutron Binding Energy

The first full experiment that was attempted was measuring the binding energy of neutrons by observing the γ -ray emitted when the neutrons are captured by protons in the water surrounding a neutron source via reaction 2.4,



The exact energy would be found from the peak from this photon observed in the spectrum.

In order to observe the spectrum as clearly as possible, the experiment was run overnight collecting data for approximately 25 hours. This would give us a much clear spectrum and so lower errors on the final data. The spectrum that was recorded is shown in figure 5, along with a description of each of the major peaks that are visible. This plot is shown on a single $\log y$ axis so that the peaks are more easily discernible.

The major peaks, from figure 5 are as follows:

1. 511 keV peak from the electron rest mass.
2. Unknown.
3. Unknown.
4. Double escape peak?.

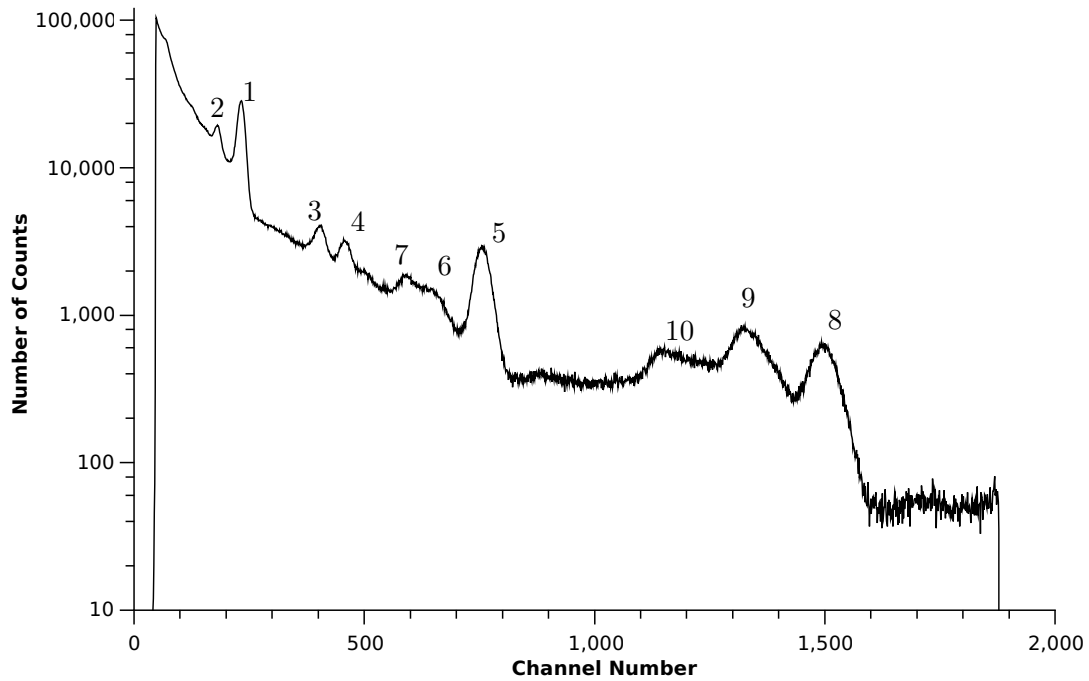


Figure 5: This is The spectrum from the neutron tank observing a radioactive neutron emitter for a prolonged period.

5. Neutron binding energy.

6. Unknown

7. Single escape peak from the neutron binding energy.

Double escape peak from the neutron binding energy.

8. Reaction peak from either oxygen ($^8\text{O}_2$) or carbon (^{12}C) from the water in the tank. Further investigation would be needed to determine the exact cause of this peak since it is so weak in comparison to the others.

9. Escape peak from the carbon/oxygen peak

10. Double escape peak from the carbon/oxygen peak

3 Neutron Attenuation In Water

3.1 Boron Tri-Fluoride Detector

The BF_3 detector contains a boron trifluoride gas at 0.5 to 1 atmospheres which acts as both a target for slow neutron conversion into secondary particles, as well as a proportional gas. A large potential difference is applied across the gas, the cathode held at ground is the outer tube of the detector and the anode, at the high voltage, is a single wire that passes through the centre of the detector, see figure 6. When a neutron is incident in the reactor, one of the reaction below takes place.



The products of these two equations differ only by the energy of the lithium which, in some cases, is left in an excited metastable state with a slightly different energy. The relative probabilities of

each reaction is 96% for equation 3.1 and 4% for equation 3.2. The resulting ions that are created are attracted to the poles of the detector and so provide a voltage through it. This provides a signal which is sent to the processing computer. It should be noted that this detector cannot give any information of the energy of the neutron since the formation of the ions is independent of energy over a certain level. Instead, this detector is used simply as a counter to obtain a value for the number of neutrons hitting it in a specified time period.

The electric field close to the wire increases quickly ($E \propto \frac{1}{r^2}$) and so a Townsend avalanche forms near to the wire when a reaction occurs. This involves the electrons being accelerated greatly by the electric field which causes the gas to conduct through the avalanche producing the signal.

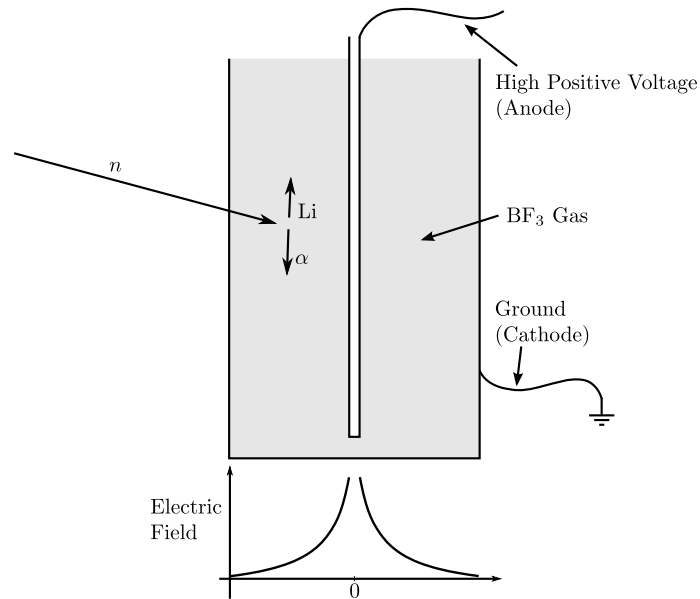


Figure 6: The BF_3 detector uses a high charge to create a flow of ions to measure the incoming neutron. It cannot measure the energy, but simply gives a count for the number of neutrons incident in the detector.

3.1.1 Wall Effect

The BF_3 detector shall be used for several purposes later, so it is worth discussing some of the features of the data that is produced by this detector. Figure 7 shows a typical spectrum from the BF_3 detector. This shows a main peak, a smaller, higher energy secondary peak and an example of the wall effect.

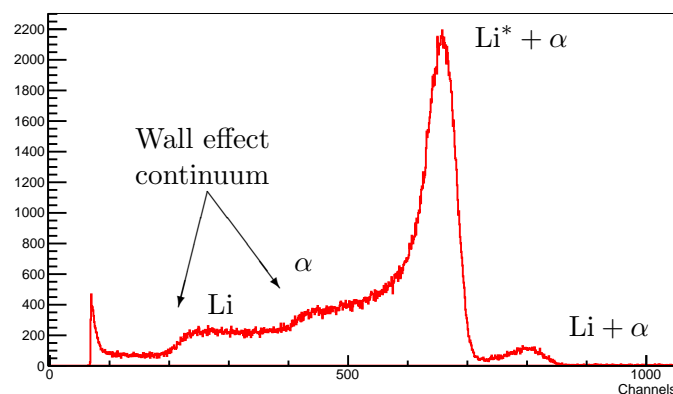


Figure 7: An exemplar spectrum from the BF_3 detector showing the wall effect when one of the products from equation 3.1 escapes the detector.

The two steps in the count number at lower channel numbers is known as the wall effect. These are due to the detector being finite and so there is a chance that one of the products from the reaction escapes from the detector. Due to the relative masses of the particles created, α and Li, the α particle takes most of the energy from the original neutron and so when it is lost from the detector, the recoded energy is lower.

Other features of the spectrum are the abrupt cut-off at the low channel numbers. This is to reduce the background counts that are present in this region from obscuring the rest of the data.